

CLIAwaived,™ Inc. INSTANT DRUG TEST CUP II



DETECTS UP TO 14 DRUGS

- Amphetamine (AMP)**
- Barbiturates (BAR)**
- Buprenorphine (BUP)**
- Benzodiazepines (BZO)**
- Cocaine (COC)**
- Ecstasy (MDMA)**
- Methadone (MTD)**
- Methadone Metabolite (EDDP)**
- Methamphetamine (MET)**
- Marijuana (THC)**
- Opiates 300 (OPI 300)**
- Oxycodone (OXY)**
- Phencyclidine (PCP)**
- Tricyclic Antidepressants (TCA)**


 12 panel includes 3 adulterants (CR, OX, pH)
 14 panel includes 6 adulterants (OX/pH/SG), (GL/NI/CR)



- FDA-cleared and CLIA-waived
- Non-negative results in five minutes
- Donor-friendly cup opening
- Fast, dark test lines
- Optional tests for adulterants
- Remarkable pricing and quality

PANELS	ADULT	AMP	BAR	BUP	BZO	COC	EDDP	MDMA	MET	MTD	OPI	OXY	PCP	TCA	THC	PACKAGING: 25 CUPS
5	■					■			■		■				■	CLIA-IDTC-5-60527-A
7					■	■			■	■	■				■	CLIA-IDTC-7-60717
12		■	■	■	■	■	■	■	■	■	■	■	■	■		CLIA-IDTC-12-NTHC
12	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	CLIA-IDTC-12-A
12 + BUP	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	CLIA-IDTC-12-BUPA
14	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	CLIA-IDTC-14-BUPA

Call your sales representative now at 1-888-882-7739



1. Remove the cup from the foil pouch.



2. Write donor ID and date on the label.



3. Collect specimen in the cup.



4. At 5 minutes, peel the label and read results. Do not read after 60 minutes.

INTERPRETATION OF RESULT

Each strip contains one drug test. C region shows validity of a test result. T region shows result for this test.

For C region:

The appearance of a line indicates a valid result.


No line means an Invalid result.

If a test strip does not have a line in the C region, test results are Invalid for this strip.

For T region:

The appearance of a line indicates a Negative result.

Note: Any test line, even a very faint test line, is considered a negative result.

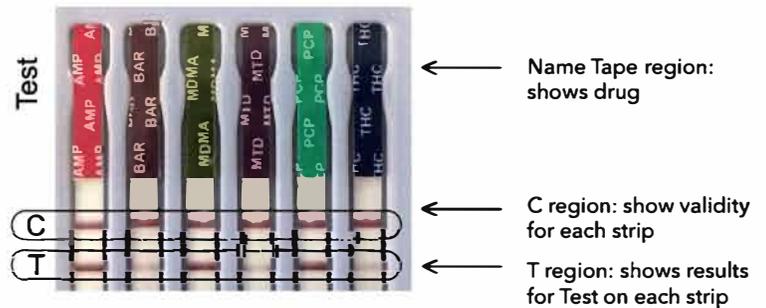

No line indicates a Preliminary Positive result.

Note: Any urine with preliminary positive results should be sent to a laboratory for confirmation.



Example #1 (Negative Result)

A red or pink line next to the "T" (drug test line) under the drug name indicates a negative result for that drug. Certain lines may appear lighter or thinner than other lines.

Therefore, it is Negative for all tests.



Example #2 (Preliminary Positive Result)

If NO red or pink line appears next to the "T" under the drug name, the sample may contain that drug.

Therefore, it is Preliminary Positive for AMP and THC.

All other tests are Negative.



Example #3 (Invalid Result)

A colored line should always appear next to the letter "C" on every test strip. If no control line appears on any of test strips, the result is invalid.

Therefore it is Invalid for both MTD and TCA tests. All other tests are Negative.