

Bayer's A1CNow+[®] Multi-Test A1C System

Distributed by

CLIAwaived[™] Inc.



For more product information visit
www.A1CNow.com

*"Study results with healthcare professionals showed that the accuracy of A1CNow+with fingerstick samples was, on average, 99%. This means that, on average, a true 7.0% A1C could read approximately 6.9% A1C. An individual A1CNow+result may differ by as much as -1.0% A1C to +0.8%A1C from the true result. This represents the 95% confidence limits of a Bland-Altman plot."

Bayer's A1CNow+[®] Multi-Test A1C System provides you with a fast and easy way of obtaining accurate A1C results in your office. This innovative technology enables you to communicate face-to-face with your patients about their diabetes control. The A1CNow+ Multi-Test A1C System provides quantitative measurement of the percent of glycated hemoglobin (A1C) levels in capillary (fingerstick) or venous whole blood samples. *(This test is for professional use only to monitor glycemic control in people with diabetes.)*

BENEFITS:

- Provides results in just 5 minutes
- Enables A1C testing in every exam room
- No maintenance or capital investment
- No refrigeration necessary if used within 4 months
- NGSP certified - National Glycohemoglobin Standardization Program (NGSP)
- 99% lab accuracy*
- FDA-Cleared / CLIA Waived

Ordering Information:

Product Description	Catalog #
Bayer A1CNow+ [®] Multi-Test A1C System (10 cartridges)	BAYER-3024
Bayer A1CNow+ [®] Multi-Test A1C System (20 cartridges)	BAYER-3021

Reimbursement Information:

CPT Code #	National Limit Amount
83037QW Hemoglobin; glycosylated (A1C)	\$14.17


DISCLAIMER – The CPT information contained in this document is provided as representative examples of reimbursement in this category. It is intended to assist providers in accurately obtaining reimbursement for health care services. It is not intended to increase or maximize reimbursement by any payor. Providers should consult their payor organizations with regard to local reimbursement policies. The information provided in this document is for information purposes only and represents no statement, promise or guarantee by CLIAwaived[™] Inc. All CPT codes are supplied for information purposes only and represent no statement, promise or guarantee by CLIAwaived[™] Inc. that these codes will be appropriate or that reimbursement will be made. CPT codes and descriptions are copyrights of the American Medical Associations. CPT does not include fee schedules, relative values or related ratings. The source for this information is the Center for Medicare and Medicaid Services. The content provided by the Center for Medicare and Medicaid Services is updated frequently. It is the responsibility of the service provider to confirm the appropriate coding required by their local Medicare carriers, fiscal intermediaries and commercial payors. All product information and prices are subject to change without notice at any time. Please call our office before ordering to verify current prices.

Bayer's A1CNow+® Multi-Test A1C System Test Instructions


Important Steps When Using The A1CNow+




1 Match lot #'s

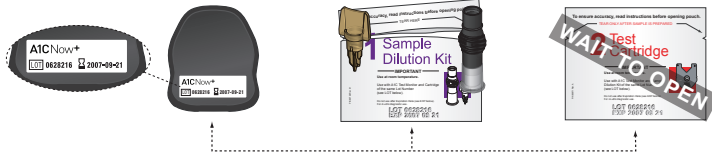


2 Collect the right amount of blood

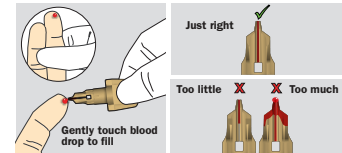

3 Insert Blood Collector all the way

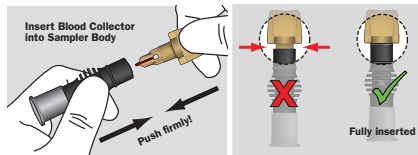

4 Shake vigorously 6-8 times

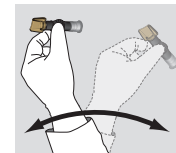

5 Insert Cartridge and Match Codes



6 Deploy Sampler and remove quickly



7 Results in 5 minutes

