


CLIAwaived™ Inc. Rapid Drug Test Cup

NEW!



The CLIAwaived Inc. Rapid Drug Test Cup is one of the simplest, easiest-to-use drugs of abuse test in the market. It provides fast, and accurate results in minutes.

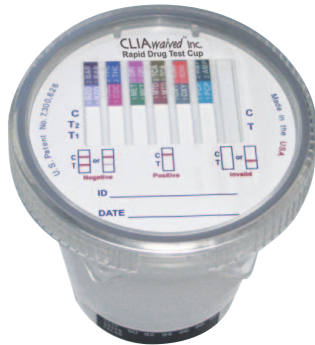

Just turn and rotate the cup to activate.

The all new CLIAwaived, Inc. Rapid Drug Test Cup (RDTC) is a CLIA Waived 12-panel urine drug testing cup. The CLIAwaived, Inc. RDTC offers instant results at point-of-care. Just turn and rotate the RDTC cup once or twice to activate the test. The cup does all the work, in only 5 minutes! Best of all, it's cost effective!

BENEFITS:

- Simple procedure
- Easy to interpret color bands
- Built-in test control
- Room temperature storage
- 14-18 month shelf life
- Integrated adulterants
- FDA-Cleared / CLIA Waived for up to 12 drugs

Ordering Information:

Product Description	Catalog #
CLIAwaived™ Inc. Rapid Drug Test Cup (THC, COC, OPI, AMP, mAMP, PCP, BAR, BZD, MTD, MDMA, TCA and OXY) Tests for common adulterant testing including: Creatinine, Nitrite, pH, Bleach, and SG.	CLIA-12 RDTC

Reimbursement Information:


CPT Code #	National Limit Amount
Effective 1/1/2011 the Healthcare Common Procedure Coding System (HCPCS) included the following new code: G0434QW – Drug screen, other than chromatographic; any number of drug classes, by CLIA waived test or moderate complexity test, per patient encounter.	\$20.83

DISCLAIMER – The CPT information contained in this document is provided as representative examples of reimbursement in this category. It is intended to assist providers in accurately obtaining reimbursement for health care services. It is not intended to increase or maximize reimbursement by any payor. Providers should consult their payor organizations with regard to local reimbursement policies. The information provided in this document is for information purposes only and represents no statement, promise or guarantee by CLIAwaived™ Inc. All CPT codes are supplied for information purposes only and represent no statement, promise or guarantee by CLIAwaived™ Inc. that these codes will be appropriate or that reimbursement will be made. CPT codes and descriptions are copyrights of the American Medical Associations. CPT does not include fee schedules, relative values or related ratings. The source for this information is the Center for Medicare and Medicaid Services. The content provided by the Center for Medicare and Medicaid Services is updated frequently. It is the responsibility of the service provider to confirm the appropriate coding required by their local Medicare carriers, fiscal intermediaries and commercial payors. All product information and prices are subject to change without notice at any time. Please call our office before ordering to verify current prices.

CLIAwaived™ Inc. Rapid Drug Test Cup Instructions

TEST PROCEDURE:

1. Urinate into the cup, filling the cup with urine above the minimum level line.
2. Secure the test lid to the specimen cup.
3. Tip the cup and rotate to activate the test.
4. Wait 5 minutes.
5. Read the results.
5. Refer to product insert for detailed instructions.

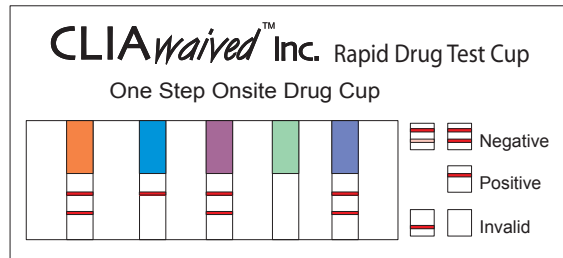

- **ACCURACY:** FDA 510(k) registered detection strips have the same high degree of accuracy when compared to SAMHSA certified laboratory screening methods.
- **SELF-CONTAINED:** Each cup includes a temperature strip and identification label. Each cup is sealed in its own bag and has an 14-16 month shelf life.
- **NO LEAKAGE:** The cup does not leak so a positive specimen can be shipped in the same container to a laboratory for GC/MS confirmation.
- **COST:** In addition to providing immediate results, the CLIAwaived Inc. Rapid Drug Test Cup costs less than half the price of laboratory testing.

ADULTERATION COLOR COMPARISON CHART

TEST AND READING TIME	ABNORMAL (LOW)	NORMAL	ABNORMAL (HIGH)
Creatinine (Cr) 45 seconds	<10 10	20 50 100 200 mg/dl	
Nitrite (Ni) 45 seconds		0 0.1-0.2 0.5-5.0	> 15 mg/dl
pH immediate	2 3	4 5 7 9	≥ 10
Bleach (Bl) 45 seconds		Negative	Positive
Specific Gravity (S.G.) 60 seconds	1.000	1.005 1.015 1.025	≥ 1.030

P/N: 88003-Rev4

Rev. 1.03.11