



Clearview®



Clearly different.  
**Clearly better.**

Women's Health | HIV / STDs | Respiratory | Cardiovascular | Enterics & GI | General Health

## Clearview® hCG

### Clearview hCG Testing

One-step pregnancy tests for early diagnosis.

### Clearly different

Clearview's patented innovative technology provides highly sensitive and specific tests for pregnancy.

- **Rapid Response:** Single-step; positive results in as little as 1 minute.
- **Simplicity:** CLIA waived (urine), easy-to-use test.
- **Reliability:** Highly sensitive assay with built-in controls ensures accuracy.
- **Flexibility:** Use with multiple sample types — urine or serum.

### Clearly better

Clearview provides better patient care with early diagnosis for proper planning, reducing potential risks to mother and fetus.

- Earlier detection allows for improved care.
- Rapid results help prevent infection or contraindicated treatment.



Point of Care Diagnostics

**CLIAwaived.com**™

CLIAwaived.com - 4332 Corte de la Fonda - San Diego, CA 92130

Toll Free : (888) 882-7739 - Phone (858) 481-5031

FAX: (801) 720-7568 - www.cliawaived.com



**Clearview hCG Combo II**  
 40 Cassettes  
 Order No. 92215  
 Urine - CLIA Waived  
 Serum - Non-waived

**Clearview 25 hCG Combo**  
 40 Cassettes  
 Order No. 92225  
 Urine - CLIA Waived  
 Serum - Non-waived

**Clearview hCG Urine**  
 50 Dipsticks  
 Order No. 92211  
 CLIA Waived

**Clearview hCG Urine**  
 40 Cassettes  
 Order No. 92217  
 CLIA Waived

**hCG Controls**  
 1 Positive, 1 Negative  
 Order No. 92218

**Clearview® hCG Combo Tests Comparison Grid**

Test Name	Clearview hCG Combo II & 25 hCG Combo	Quidel QuickVue® One-Step hCG Combo	Beckman Coulter ICON® 25 hCG	Genzyme OSOM® hCG Combo Test
Format	40 Cassettes	50 Cassettes	25 Cassettes	25 Cassettes
CLIA Status: Sample Types	Waived: Urine Non-waived: Serum	Waived: Urine Non-waived: Serum	Waived: Urine Non-waived: Serum	Waived: Urine Non-waived: Serum
Sensitivity	<b>hCG Combo II</b> 20 mIU/mL Urine 10 mIU/mL Serum	25 mIU/mL	25 mIU/mL	20 mIU/mL Urine 10 mIU/mL Serum
	<b>25 hCG Combo</b> 25 mIU/mL			
Specificity	>99%	>99%	>99%	>99%
Test Steps	1	1	1	1
Minutes to Result	(+) = 1* (-) Urine = 3 (-) Serum = 5	Urine = 3 Serum = 5	Urine = 3 Serum = 5	Urine = 3 Serum = 5
Shelf Life	24 months	18 months	24 months	24 months
Storage	Room Temperature	Room Temperature	Room Temperature	Room Temperature

**Clearview hCG Urine Dipstick Procedure**




Dip to line,  
5 seconds



Read in  
3 minutes

Interpret results

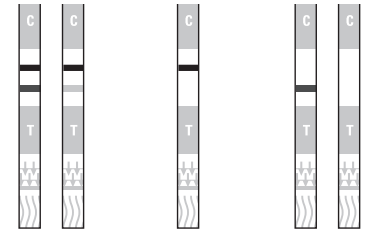

**Positive**


**Negative**

**Invalid**

**Clearview hCG Combo II / 25 Combo / Urine Cassette Procedure**



Draw sample

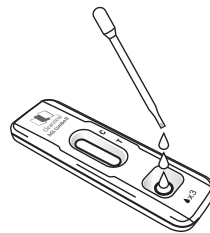

Add 3 drops  
to sample well



Read in  
3 minutes - Urine  
5 minutes - Serum



**Positive**

**Negative**

**Invalid**

Interpret results

visit

[www.invernessmedicalpd.com](http://www.invernessmedicalpd.com)  
 for more information



**CLIAwaived.com™**

CLIAwaived.com - 4332 Corte de la Fonda - San Diego, CA 92130

Toll Free : (888) 882-7739 - Phone (858) 481-5031

FAX: (801) 720-7568 - [www.cliawaived.com](http://www.cliawaived.com)