



Clearview®



Clearly different.
Clearly better.

Women's Health | HIV / STDs | Respiratory | Cardiovascular | Enterics & GI | General Health

Clearview® hCG

Clearview hCG Testing

One-step pregnancy tests for early diagnosis.

Clearly different

Clearview's patented innovative technology provides highly sensitive and specific tests for pregnancy.

- **Rapid Response:** Single-step; positive results in as little as 1 minute.
- **Simplicity:** CLIA waived (urine), easy-to-use test.
- **Reliability:** Highly sensitive assay with built-in controls ensures accuracy.
- **Flexibility:** Use with multiple sample types — urine or serum.

Clearly better

Clearview provides better patient care with early diagnosis for proper planning, reducing potential risks to mother and fetus.

- Earlier detection allows for improved care.
- Rapid results help prevent infection or contraindicated treatment.



Point of Care Diagnostics

CLIAwaived.com™

CLIAwaived.com - 4332 Corte de la Fonda - San Diego, CA 92130
Toll Free : (888) 882-7739 - Phone (858) 481-5031
FAX: (801) 720-7568 - www.cliawaived.com



Clearview hCG Combo II
40 Cassettes
Order No. 92215
Urine - CLIA Waived
Serum - Non-waived

Clearview 25 hCG Combo
40 Cassettes
Order No. 92225
Urine - CLIA Waived
Serum - Non-waived

Clearview hCG Urine
50 Dipsticks
Order No. 92211
CLIA Waived

Clearview hCG Urine
40 Cassettes
Order No. 92217
CLIA Waived

hCG Controls
1 Positive, 1 Negative
Order No. 92218

Clearview® hCG Combo Tests Comparison Grid

Test Name	Clearview hCG Combo II & 25 hCG Combo	Quidel QuickVue® One-Step hCG Combo	Beckman Coulter ICON® 25 hCG	Genzyme OSOM® hCG Combo Test
Format	40 Cassettes	50 Cassettes	25 Cassettes	25 Cassettes
CLIA Status: Sample Types	Waived: Urine Non-waived: Serum	Waived: Urine Non-waived: Serum	Waived: Urine Non-waived: Serum	Waived: Urine Non-waived: Serum
Sensitivity	hCG Combo II 20 mIU/mL Urine 10 mIU/mL Serum	25 mIU/mL	25 mIU/mL	20 mIU/mL Urine 10 mIU/mL Serum
	25 hCG Combo 25 mIU/mL			
Specificity	>99%	>99%	>99%	>99%
Test Steps	1	1	1	1
Minutes to Result	(+) = 1* (-) Urine = 3 (-) Serum = 5	Urine = 3 Serum = 5	Urine = 3 Serum = 5	Urine = 3 Serum = 5
Shelf Life	24 months	18 months	24 months	24 months
Storage	Room Temperature	Room Temperature	Room Temperature	Room Temperature

Clearview hCG Urine Dipstick Procedure




Dip to line,
5 seconds



Read in
3 minutes

Interpret results

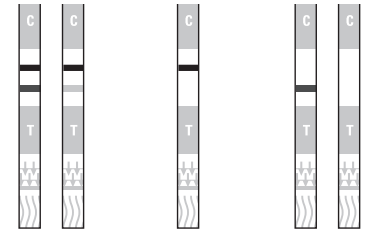

Positive


Negative

Invalid

Clearview hCG Combo II / 25 Combo / Urine Cassette Procedure



Draw sample

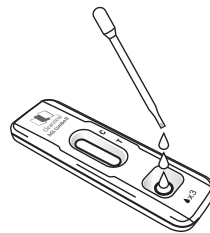

Add 3 drops
to sample well



Read in
3 minutes - Urine
5 minutes - Serum



Positive

Negative

Invalid

Interpret results

visit

www.invernessmedicalpd.com
for more information



CLIAwaived.com™
 CLIAwaived.com - 4332 Corte de la Fonda - San Diego, CA 92130
 Toll Free : (888) 882-7739 - Phone (858) 481-5031
 FAX: (801) 720-7568 - www.cliawaived.com